ANLAN

Always for you

www.anlan.com

Designed by ANLAN

Made in China

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ANLAN Guasha Beauty Device

USER MANUAL

※ Thank you for purchasing the ANLAN Guasha Beauty Device. Before using this product, please read this user manual carefully for future reference. In addition, this manual also serves as a warranty, so please keep it properly.

Warranty Certificate

Product Name	ANLAN Guasha Beauty Device
Model Type	DYM-078
Warranty Period	1 Year Warranty
Customer Information	Name:
	Phone Number:
	Address:
	Email:
Place of Purchase	

Safety Precautions

Please read the document content carefully to make sure you understand all information.

Please be sure to comply with this precautions to prevent users and others from being injured or property loss.

	Indicates that incorrect handling may result in death or serious injury to the user.
	Indicates that incorrect operation may result in minor injury to the user or property damage.

		Danger
	Do not use the following medical electrical equipment with this device together:	
	• Implantable medical electronic devices such as pacemakers • Medical electronic devices for life-threatening items such as heart-lung machines • Wearable medical electronic devices such as electrocardiographs	
	Do not use this product if you are in such situation :	
	• Receiving treatment from a doctor. • During pregnancy (can use at least 6 months after childbirth, please consult a doctor before use.) • Menstrual period, lactation Minors, and people who cannot express their thoughts clearly.	

	• Minors, and people who cannot express their thoughts clearly. • People with cancer, epilepsy or other functional disabilities • People with blood diseases, infectious diseases, hypertension, malignant tumors, visceral diseases, acute diseases, tuberculosis, and fever. • People with allergies, sensitive, fragile skin, atopic dermatitis, and skin problems. • People with skin diseases or heat allergies • People who have developed lesions due to physical stimulation, etc. • Do not use this product near equipment that emits short waves or microwaves. • People who are undergoing dental treatment or fillings may experience tooth discomfort due to EMS irritation. • Please lower the level and adjust the area you want to use. ※ The above situations may cause accidents or malfunctions.
---	--


	 Note
	Do not use it in the following cases:
	● If you are extremely fatigued or in poor physical condition. ● While taking steroid hormones, immunosuppression, or drugs that cause drowsiness. ● Areas where the skin has acne, cuts, bulging moles, viral warts, or areas where the skin is suppurated and inflamed. ● Areas where the skin is irritated and abnormal due to excessive sunburn. ● Area where hair is being removed or removed on the day of use. ● If you have undergone facial surgery (plastic surgery, cosmetic surgery, etc.). ※ The above may cause accidents or trouble.

Warning

※ If you are unsure, please consult your doctor beforehand to ensure that the beauty device can be used, and if you have any questions, please contact service. For detailed usage precautions, the contents listed in the instruction manual take precedence.

Product Introduction

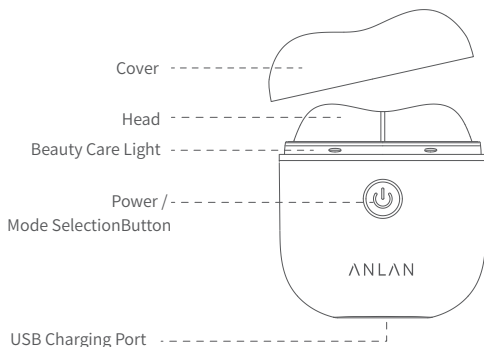※ No operation in each mode, automatically shuts off after 10 minutes.

Pack List

Main Device × 1

Instruction Manual × 1

Type-C Cable × 1

Specifications

Product Name	ANLAN Guasha Beauty Device
Model Type	DYM-078
Model Number	DYM-078-PK
Size	About 80 × 75 × 26 mm
Weight	About 84g
Material	ABS + Zinc Alloy
Temperature	Approx. 38 ± 0.7°C – 46 ± 0.7°C
Input	5V=1A
Battery Capacity:	400mAh
Auto Power-Off Time	About 10 minutes
Charging Time	About 70 minutes
Continuous Operating Time	About 45–80 minutes
LED Beauty Light	7 Color Light Therapy (Red, Cyan-Green, Yellow, Green, Blue, Purple, White)

About Modes & Operation

Modes	Function
Mode 1	Red Light+Heating (38 ± 2°C)
Mode 2	Cyan-Green Light+EMS (Low)
Mode 3	Yellow Light+Vibration (Low)
Mode 4	Green Light+EMS (Medium)+Vibration (Medium)
Mode 5	Blue Light+EMS (Low)+Heating (42 ± 2°C)
Mode 6	Purple Light+EMS (Medium)+Vibration (High)
Mode 7	White Light+EMS (High)+Heating (46 ± 2°C)

Operating Instructions

● Power On

Press and hold the Power / Mode Selection Button () to turn the device ON.

※ Auto Power-Off Function: If no button operation is performed for 10 minutes, the device will automatically power off.

● Adjust the Mode

Press the Power / Mode Selection Button () briefly to switch modes in the following order:

Mode 1 → Mode 2 → Mode 3 → Mode 4 → Mode 5 → Mode 6 → Mode 7.

※ While the device is powered on, each short press of the Power / Mode Selection Button () will switch to the next mode, and the corresponding beauty light will illuminate.

● Power Off

After use, press and hold the Power / Mode Selection Button () to turn the device OFF.

How to Use

Note:


- Be sure to clean the skin in the massage area before use
- Using micro-current mode requires keeping the skin moisturized to a certain degree
- If your skin becomes dry during use, please replenish and apply moisturizing products appropriately in order for the EMS micro-current to work better.
- You can follow the recommended usage steps, or you can combine modes at will or use a mode alone, depending on your preference.

For Face Care

1. Follow your daily routine and thoroughly clean your face.
2. Apply skin care products such as toner, essence, lotion on face.

3. Use Mode 1 to lift face upward to promote absorption.

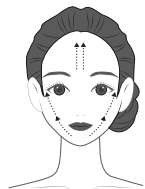
4. After the essence is absorbed, use Mode 2 to lock in nutrients.
5. Use Mode 3 to massage upwards to improve edema, tighten skin and reduce wrinkles

For Body Care


1. Take a shower or wash body skin before use
2. Apply body lotion, gel or essential oil on back, arm or leg.


3. Choose the modes you like to massage skin

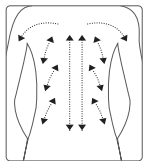Back

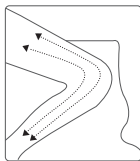Leg

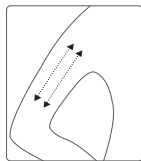Arm

Charging Instructions

Please fully charge the machine before using it for the first time.

- If you have not used this product for a long time, please fully charge it before use

In order to protect the built-in battery, it needs to be charged once a month even if it is not used for a long time.

- When the battery is low, connect the machine to the charging cable and charger, and plug it into a household power socket;
- During charging, The power button LED will blink while charging. When fully charged, the power button LED will stop blinking.
- Too low or too high temperature will shorten the battery life, please charge at room temperature of 15°C to 30°C. The main device may heat up slightly while charging, but this is not a malfunction.

- ※ The device cannot be used while charging.
- ※ Before charging, please make sure the charging port is clean and dry.
- ※ Please be sure to use the correct voltage adapter mentioned in this manual for charging.

About Cleaning And Storage

- Please turn off the machine in time after use. If any lotion, essence, or gel remains on the machine, please clean it in time.
- Be sure to turn off the power when cleaning to avoid electric shock, short circuit or fire.
- This product is not waterproof, please do not wash it with water to avoid machine damage and functional failure.
- Wipe away dirt with a soft cloth dampened in water and wrung out. If the dirt does not come off, wipe it gently with a cloth dampened with dilute neutral detergent.
- Do not use alkaline detergents, alcohol, thinner, gasoline, nail polish remover, etc.
- The use of oily substances, bleach, etc. is prohibited.
- Do not disassemble, repair or modify it by yourself.
- The warranty excludes injuries, malfunctions, damages caused by incorrect use, disassembly, modification, repair, etc., and we are not responsible for them once the above precautions and warnings are not followed.

Frequently Asked Questions (FAQ)

Q: The cap is difficult to remove.

A: Firmly hold the cap with your hand and pull it downward to remove.

Q: I cannot feel the EMS.

A: EMS transmission may be poor if the skin is too dry, which can reduce effectiveness or cause skin issues. If EMS feels weak, apply a generous amount of gel or lotion to ensure the entire electrode head is in contact with the skin.

※ Creams or oils containing fats may interfere with EMS and are not recommended.

Q: Can I use other gels or cosmetics?

A: We recommend the dedicated gel for optimal EMS sensation, but toner or serum can also be used.

Q: Can I use it in the bath?

A: The device is not waterproof. Do not use it in the bathroom or wet areas.

Q: Can I use it continuously?

A: We recommend using the device once per day for about 5 minutes. If using consecutively, allow an interval of at least 5 hours between sessions.

Q: I feel a tingling sensation.

A: The device emits a microcurrent to support skin care, which may cause mild stimulation. For first-time use, start at Level 1. Apply gel to ensure full contact between the electrode head and the skin.

JP

JP

EN

EN

※ Stop using immediately if acne develops, the skin becomes sensitive, or you feel strong discomfort. Consult a dermatologist if necessary.

Q: The device cannot be charged.

A: Check the charging port for debris and clean if necessary. Try using another cable and a 5V/1A charger. If the issue persists, contact customer support.

※ Apple charging adapters and Type-C to Type-C cables with MFI certification cannot be used to charge this device.

Product Disposal

JP

EN

This product contains a lithium-ion battery.

Lithium-ion batteries are valuable and recyclable resources. When disposing of the product, please take it to a cooperating store with a recycling box or to a designated collection point.