
ANLAN

Always for you

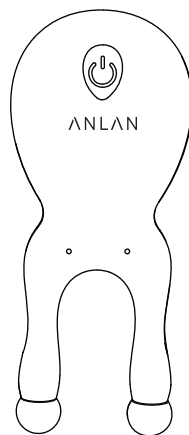www.anlan.com

Designed by ANLAN

Made in China

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
ANLAN

ANLAN Light Frog Facial Massager User Manual

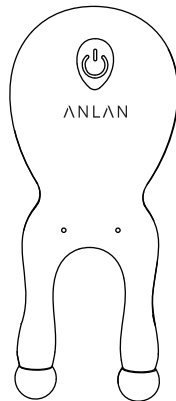Thank you for purchasing the ANLAN Light Frog Facial Massager. Before using this product, please read this user manual carefully for future reference. In addition, this manual also serves as a warranty, so please keep it properly.

Warranty Certificate

Product Name	ANLAN Light Frog Facial Massager
Model	DYM-070
Warranty Period	1 Year Warranty
Customer Information	Name:
	Phone Number:
	Address:
	Email:
Place of Purchase	

Safety Precautions

Important safety instructions

When using electrical appliances, especially when children are present, be careful and follow the following safety instructions:

- Please read the instructions carefully before use
- The appliance should be kept away from water sources

Safety precautions

WARNING

Please be careful with plastic bags used on packaging or products, and keep them away from children. Do not use in cribs, beds, strollers or children's play areas, Because the thin layer of the plastic bag may block the nose and mouth, hindering breathing. Please note that this plastic bag is not a toy.

WARNING

Do not use this appliance near a bathtub, shower, wash basin or other container filled with water.


- Please use the charging cable randomly provided with this product or a charging cable that meets the product specifications, and correctly connect the charger plug to the household socket to avoid malfunctions, fires, or fires caused by heat.

- Due to the built-in rechargeable battery of this product, please do not place it in a fire or use it in high temperature or high heat environments, as this may cause malfunctions, fires, or fires due to heat generation.
 - If the product heats abnormally, deforms, damages, or has an odor, do not use it again to avoid the risk of electric shock and fire.
 - This product is a personal massage product and cannot replace drugs or medical equipment.
 - This product is not suitable for use by individuals with physical, sensory, or mental problems, or those lacking experience and common sense (including children), unless supervised and guided by a guardian, and responsible for the safety of the user.
 - Supervise children to ensure they do not play with this tool.
 - Keep packaging materials (such as plastic bags) out of reach of children.
- If the machine is damaged or the product functions abnormally, in order to avoid danger, it must be repaired by professionals from the manufacturer, its maintenance department, or similar departments.
- This product is not waterproof, so it is important to prevent water from splashing onto the surface of the product during use. Please do not use this product while taking a shower.
 - When cleaning and maintaining, please turn off the machine and perform it while it is turned off.
 - Do not store or use this product in the bathroom or other damp rooms.
 - Do not place this product near heat sources or flammable areas, as it is prone to malfunction and danger, and may cause deformation, discoloration, or damage.

- If you feel any abnormality with this device during use, please stop using it immediately and contact the factory's after-sales service.
- If you experience any physical abnormalities or pain during the use of this device, please stop using it immediately and consult a healthcare professional. People with dermatitis, skin damage, eczema, swelling, and other conditions are not recommended to use it.
- This product belongs to tool care equipment, not medical devices, and does not have therapeutic effects. If you need treatment or medication, please consult a professional doctor or relevant personnel.

Introduction to Product Functional Parts

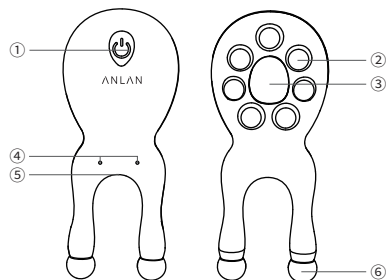1. Power ON/OFF, Mode Switching
2. Massage Ball
3. Beauty Lamp
4. Pilot Lamp
5. Charging Port
6. Massage Ball

Specifications

Product Name	ANLAN Light Frog Facial Massager
Model Type	DYM-070
Model Number	DYM-070-PK
Battery Life	About 100 Minutes
Battery Capacity	250mAh
Net Weight	Appr. 74g
Input Voltage And Current	5V --- 1A
Charging Time	About 120 minutes
Dimensions	Appr. 113 x 47 x 26mm
Material	ABS, PC, Zinc Alloy
Flash Energy/J	7modes (Red light, Cyan light + EMS, Yellow light + vibration, Green light+EMS, Blue light + EMS + vibration, Purple light + EMS + vibration, White light + EMS + vibration)

Packing List

Host * 1

Type-C
Charging Cable * 1

Instruction
Manual * 1

Product Usage


1. Wash your face with a detergent and wipe your face with a towel.

2. Apply tonic/balm/serum etc.

3. Use roller on face, neck, shoulders.

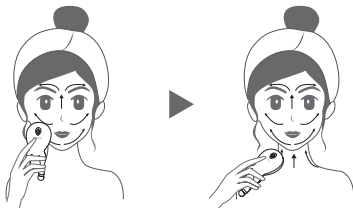
4. Follow the instructions to raise the bloating line and shape the contour.

- ① From chin to ear along jaw line
- ② From the corner of the mouth to the cheekbones
- ③ From nose to forehead
- ④ From the middle of the face to the outside
- ⑤ From the inside of the forehead to the outside

Maintenance

Before cleaning, please turn off the machine.

- Protect this device from moisture.
- Do not use sharp or corrosive or volatile diluents for cleaning. You can use a damp soft cloth or cotton pad to gently wipe.
- Do not rinse this device directly with water to avoid damage caused by water entering the machine.

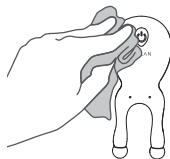Warranty

1. We provide a one-year warranty for our products, which starts from the date of purchase.
2. During the warranty period, if the product has defects due to materials or production processes, we will replace the parts or the entire product as appropriate.
3. To apply for warranty, the product and proof of purchase date (purchase invoice) must be submitted to an authorized repair station or factory after-sales service.
4. The company does not provide warranty in the following situations:
 - Damage caused by improper use, such as using the wrong power source, inserting the product into an inappropriate socket, human damage, etc.

- The aging and wear of the product caused by normal use have little impact on its functionality.
- Product problems caused by improper use for other purposes, improper operation, non professional maintenance, or the use of non company accessories
Consumables and accessories such as consumables and complimentary accessories.
- Damage caused by force majeure factors such as lightning, earthquakes, fires, floods, etc.
- After sales service does not include damages caused by the aforementioned reasons.
- This device is equipped with a lithium battery inside, so:
Before discarding the appliance, the battery must be removed from the appliance;
- When removing the battery, the appliance must be turned off and powered off;
- Batteries should be disposed of safely;
- Users are not allowed to replace batteries on their own.
- We suggest that you send the product to the official collection point and have professionals remove the battery.

— JP
—
EN