
ANLAN

Always for you

www.anlan.com

Designed by ANLAN

Made in China

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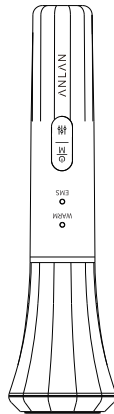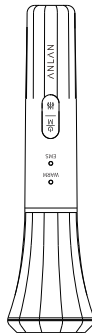
JP

EN

※ This product is not a medical device used for diagnosis or treatment, but a general household item for eyes, face and neck.

Thank you for purchasing the ANLAN EMS Eye Beauty Device. Please read this entire manual carefully before using it to ensure proper use and safety. In addition, this manual is also available as a guarantee, please safe keeping.

Warranty Certificate

Product Name	ANLAN EMS Eye Beauty Device
Model Type	AMYY31
Warranty Period	1 Year Warranty
Customer Information	Name:
	Phone Number:
	Address:
	Email:
Place of Purchase	

JP
EN

Safety Precautions

In order to prevent harm and property damage to users, please read this manual carefully and follow the manual step.


 Warning	It means inappropriate use can lead to undesirable consequences.
 Attention	It means you must follow the instructions to use the machine to avoid injury.

 Prohibition	 Warning
	Never use this device with following medical electrical equipment. • Tiny electro-magnetically damaged implant, such as pacemaker. • Medical electronics for life support such as artificial heart / artificial lung. • Portable medical electronics, such as device of electrocardiogram.

JP
EN

Prohibition JP EN	⚠ Attention
	Never use under such situation:
	● Skin with wound ● Skin area with plastic surgery ● Eye area with disease ● Inflamed skin ● Eyeball
	⚠ Warning
	People mentioned below should consult a doctor before using the device.
	● People with malignant tumors ● People with blood disorders such as bleeding disorders and hemophilia ● People with skin diseases, atopic dermatitis, facial neuropathy ● People with allergies, especially those with sensitive skin ● People who cannot feel temperature ● People who have skin problems due to excessive sun damage

Product Details

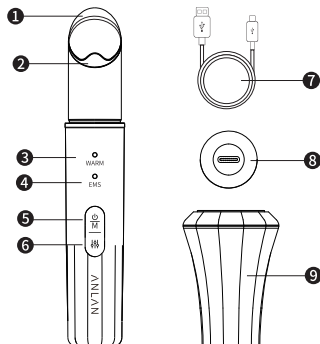- ① Massage Head
- ② Red/Blue Light
- ③ Warm Indicator
- ④ EMS Indicator
- ⑤ Power/ Mode Button

- ⑥ Level Button
- ⑦ Type-C Charging Cable
- ⑧ Charging Port
- ⑨ Dust Cover

Specification

Product Name: ANLAN EMS Eye Beauty Device
Model Type: AMYY31
Model Number: AMYY31-PK
Material: ABS
Weight: About 55g
Size: D37mm*H132mm
EMS Level: 3 Levels(indicator lamp white-green-blue)
Light: Red light (About 610nm), Blue light (About 460nm)
Temperature: About 40°C
Battery Capacity: 400mAh
Input: 5V=1A
Voltage: 3.7V
Automatic Shutdown: About 5 mins
Charging Time: About 2.5 hours
Working Time: About 80 mins
Vibration Frequency: 8000rpm

Packing contents

Eye Beauty
Device*1

User manual
*1

Type-C charging
cable*1

Precautions Before Use

1. If you wear contact lenses, be sure to take them out before use.
2. Wash your face before use to remove dirt such as makeup and sebum and apply some lotion or serum to your skin after cleansing.
3. Never use if there are any inflammation or abnormalities on your skin due to acne or pimples.
4. Do not use it continuously on the same part of the skin but move it while using. If it feels too hot, stop using it immediately.
5. Do skin test before using for the first time to make sure you are not allergic to metal.
6. When using the vibration function, never use it on your eyelids to prevent any eye problems and accidents.

Button description

Power /Mode Button:

Long press the button to turn on the power.
The default setting is red light beauty treatment / vibration / WARM mode.
If you press this button for a short time each time, the mode will switch in the order of "red light beauty treatment / vibration / WARM mode" → "blue light beauty treatment / EMS mode" → "blue light beauty treatment / EMS mode" → "blue light beauty treatment / EMS mode".
Long press the power button to turn it off.

red light beauty
treatment / vibration
/ WARM mode

blue light beauty
treatment / EMS
mode

Level Button:

PLEASE NOTE: Warm mode does not have the level, only EMS mode has.
Change to EMS mode, the default setting is level (The EMS indicator lamp is white.) Short press this button to switch (white-green-blue) "weak" → "medium" → "strong" →.

weak
(white)

medium
(green)

strong
(blue)

Instructions


Instructions for using eye care

- ① After applying eye cream or lotion, move the eyelid back and forth from the lower to the upper eyelid while drawing a circle.
- ② Move from the lower and upper eyelid to the temple area in a semicircle.

Whole face usage

- ① After washing your face, apply skin lotion, serum, etc., and then select your preferred mode.
- ② Move the massage head from the mouth and chin along the face line and lift the cheeks.

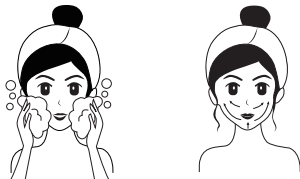Note: If you can't feel the EMS, your skin may be dry. Be sure to apply gel or lotion before using.

Charging Method

1. When using for the first time after purchase or when not in use for a long time, please fully charge the main unit.
2. Insert charging cable into the port. You can connect it to a USB port on a computer, power bank, or adapter, etc.
3. Cannot be used while charging. Unplug the cable from the main unit when not charging.
※The main unit may become slightly warm during charging, but this is not a malfunction.
4. If the battery is low, the WARM and EMS indicator (white) will blink and turn off automatically after a few seconds.

5. The WARM and EMS indicator (white) flash when charging, and it takes about 2.5 hours to fully charge. WARM and EMS light (white) light s up when charging is complete.
6. The continuous use time is about 80 minutes. Auto Power Off Function: Auto power off after 5 minutes of no operation or 5 minutes of use in each mode.

Maintenance And Storage

1. The main body is not waterproof. Do not wash with water.
2. After use, wipe off gel, lotion, serum, eye cream, etc. on the head and body with a soft cloth or tissue.
3. Do not leave the main unit in a humid place such as a bathroom. (Causes of mold and bacteria growth/Causes of failure due to condensation and rust inside the main unit)
4. Do not use alcohol, thinner, benzine, nail polish remover, etc. (Cause of faulty, cracking, discoloration of parts, etc.)
5. This product is only used as an auxiliary device. If you feel any abnormality on your skin, discontinue use and consult a specialist doctor.
6. The main body is not a toy. Keep out of reach of children.

About the Disposal

This product uses a lithium-ion battery. Lithium-ion batteries are valuable recyclable resource. When disposing of used products, remove the lithium-ion battery and take it to a cooperating store that has a recycling box or a specific collection window.

JP

EN