ANLAN

Always for you

www.anlan.com

Designed by ANLAN

Made in China

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ANLAN FACE & NECK BEAUTY DEVICE

USER MANUAL

This product is not a medical device used for diagnosis or treatment, but a general household item. Please read this user manual carefully before use. Please keep this user manual after reading for future reference. Please attach this user manual when transferring this product to others.

Warranty Certificate

Product Name	ANLAN FACE & NECK BEAUTY DEVICE
Model Type	ES-081
Warranty Period	1 Year Warranty
Customer Information	Name:
	Phone Number:
	Address:
	Email:
Place of Purchase	

SAFETY PRECAUTIONS


We do not recommend using this product for the following people:

- ※ People who have facial plastic surgery
- ※ People with metal implants under the skin
- ※ Those who are pregnant

PRECAUTIONS FOR USE

1. This product can be used to care for the skin on the face and neck.
2. It is recommended to use this product for no more than 20 minutes a day. It is not recommended for sensitive skin to use it every day.
3. Press and hold the power button for **more than 1 second** to turn on the phone, and short press to switch modes.
4. About vibration: This product uses contact vibration. **Your fingers need to hold the sensing strips on both sides of the machine** to form a closed loop with the body.
5. **About charging: The red light flashes continuously during charging, and the red light is always on after full charge. The charging time is 2 hours, and the machine cannot be turned on during charging.**
6. **The device has a safety setting of automatic power off in 5 minutes. If you need to use it again, please restart it.**
7. Please use facial cleanser to remove makeup or dirt on your face or neck first, and then use CLEAN mode to deeply cleanse the skin through the ion + export function.
8. When using EMS mode and moisturizing mode, it is recommended to use moisturizing or functional gel, lotion, essence, etc. to make EMS more obvious and more effective.
9. Please wipe the beauty head with paper. Do not rinse with water.
10. Do not place it in a humid place such as a bathroom, which may cause a short circuit.



PRODUCT INTRODUCE

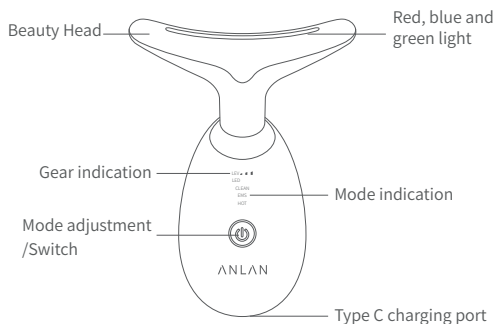SPECIFICATIONS

Product Name	ANLAN FACE & NECK BEAUTY DEVICE
Model Type	ES-081
Model Number	ES-081-PK
Material	ABS
Output voltage	3.7V
Battery capacity	500mAh
Power	2.5W
Input voltage	5V=1A
Size	89x47x120mm
Auto turn off	About 5 min
Charging time	About 2h
Package including	Face & Neck Beauty Device, Type C Charging Cable, User Manual

CHARGING

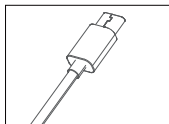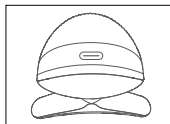1. Type C charging port, you only need to use a little effort to connect it when charging;
2. Charging response: The red light flashes continuously during charging, and the red light is always on after it is fully charged;
3. The charging time is 2h, and the device cannot be turned on during charging.

20>

FUNCTIONS

This product is a color light neck beauty instrument that can clean and brighten skin, EMS care, warm and moisturize, and integrate positive ions + negative ions + high frequency vibration + red light + blue light + green light.

Positive and negative ions

Utilize the principle of the absorption of positive and negative ions to remove the dirt in the deep pores and thoroughly clean them; through iontophoresis, the beauty ingredients are penetrated into the skin.

Vibration

Use mechanical vibration and voltage stimulation to tighten skin, increase elasticity, and increase metabolism.

Note: This face & neck beauty device uses contact vibration. When using it, when the massage head touches the skin, your hands must also hold the metal sensing strips on both sides of the product to form a closed loop to generate vibration.

Warm

Warm massage promotes absorption, so that wrinkles can be lightened and skin will become glossy.

Red light

It can effectively stimulate the skin, revitalize collagen, improve dark skin and reduce fine lines and wrinkles.

21>

Blue light

It can effectively improve skin problems, reduce acne and pimple and help to kill bacteria. Long term use can rejuvenate skin and lead to a health skin.

Green light

It can promote the microcirculation of the skin surface, dredge lymph and reduce swelling.

FEATURES

Clean mode (CLEAN)

Clean up trash from deep pores to shrink pores, eliminate acne, and solve dull skin.

EMS mode (EMS)

This mode can effectively stimulate and relax the skin and help skin care product go deeply to keep skin moist.

Moist mode (HOT)

This mode gently awakens and relaxes the skin, helping the skin restore its vitality.





HOW TO USE





Note: The device has a safety protection setting. **The product will automatically power off after 5 minutes of use. If you need to use it again, just restart it.**

Clean mode

- ※ This mode can effectively clean the dirt in the pores and brighten the skin.
- ※ Used with green light, it can promote the microcirculation of the body surface, dredge lymph and reduce swelling.
- ① Long press the power button more than 1 second to turn on with "CLEAN" mode (this mode is worked with green light)
- ② Spread cleansing water evenly on the skin of the face or neck (the machine is not waterproof, Please do not pour water directly on the machine)
- ③ **Hold the metal sensor strips on both sides of the product** (to form a cycle with the human body to achieve real ion import and export)
- ④ Move the beauty head from bottom to top of the neck or face to lift the skin.
- ⑤ Clean your neck with cleanser and water after using clean mode.

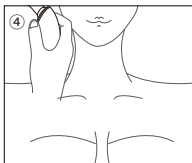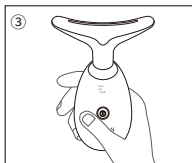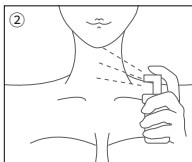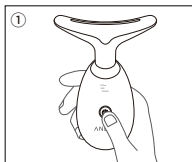- ② Spread the usual toner or lotion evenly on the skin.
- ③ **Hold the metal sensor strips on both sides of the product.**
- ④ Move the beauty head from bottom to top of the neck or face to lift the skin.

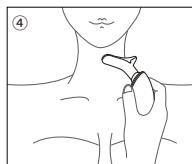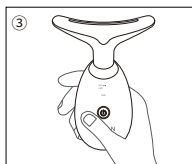
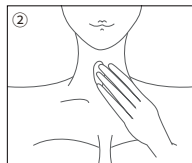
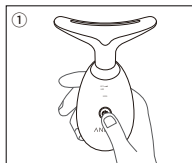
EMS mode

- ※ This mode can effectively stimulate and relax the skin and help skin care product go deeply to keep skin moist.
 - ※ Used with red light, it can effectively stimulate the skin, revitalize collagen, improve dark skin and reduce fine lines and wrinkles.
- ① Long press the power button more than 1 second to turn on, and short press the power button again to switch to the "EMS" mode (this mode is worked with red light)

Moist mode

- ※ This mode can gently awakens and relaxes the skin, helping the skin restore its vitality.
 - ※ Used with blue light, it can effectively improve skin problems, reduce acne and pimple and help to kill bacteria. Long term use can make skin rejuvenation and lead to a health skin.
- ① Long press the power button for 1 second to turn on, and short press the power button twice again to switch to the "Moist" mode (this mode is matched with blue light)

- ② Apply facial mask or neck mask.
- ③ **Hold the metal sensor strips on both sides of the product.**
- ④ Move the beauty head from bottom to top of the neck or face to lift the skin.

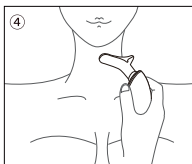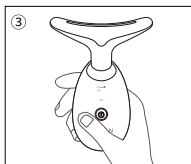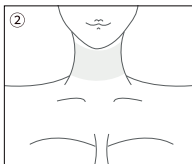HOW TO MOVE THE BEAUTY HEAD


- ① Lift from the chin to the cheek.

- ② Lift from the cheek to the corner of the eye.

- ③ Lift from the lower neck to the chin.

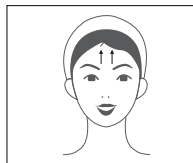- ④ Lift from the lower forehead to the upper forehead.

Q1: How to solve the problem that the product is difficult to charge?

A1: First of all, there is nothing wrong with the product charging cable. You only need to use a little force to connect it to the socket when charging.
Secondly, if you still cannot charge successfully after trying multiple methods, please contact us and we will help you solve it as soon as possible.

Q2: How to judge whether the product is fully charged? How long does it take to fully charge? How long can it be used after it is fully charged?

A2: Please charge the product before using it; the red light flashes during charging, and the red light is always on after it is fully charged. The charging time of the product is 2 hours; it can be used for 2-3 weeks after being fully charged.

Q3: Why does the product automatically shut down after a period of use?

A3: This product has a safety setting of automatic power-off. The product will automatically shut down after 5 minutes of continuous use. If you need to use it again, please restart it.

Q4: Why does the facial massager sometimes not vibrate?

A4: First of all, when the massage head touches the skin, the hands must also hold the sensor strips on both sides of the product to produce vibration. Secondly, it is necessary to keep the face moisturized. Too dry face will also cause the vibration to disappear or the vibration to be not obvious.

Q5: What skin type is suitable for? What products need to be used together?

A5: This product is suitable for all skin types. Do not use it too frequently on sensitive skin. Do not use it when the face is injured; it is recommended to apply gel, essence lotion, etc. before use for better results.

Q6: How long will it take to see the effect?

A6: First, this facial massager needs to be used regularly for more than 28 days. Second, it needs to be used with functional products such as gel and essence. Third, do not use different modes every time. We recommend that you choose the corresponding mode for the problem you need to solve first, and wait for the problem you need to solve to be improved before choosing other modes. This will make it easier to see the effect.

Q7: How to return the product?

A7: First, you can initiate a return and refund application on the order details page. After Amazon approves your application, it will give you instructions for the return.
Second, since return and refund are Amazon's system functions, we have no right to manage them. If you encounter any problems with returns, please contact Amazon customer service at the front desk.